

CAPTURE LIFE'S ADVENTURES.

Choose the Smart CPAP Choice

The 365 miniCPAP is a completely integrated and portable CPAP & heated humidifier that's ready to go wherever so you can sleep comfortably every night, everywhere.



PN 400206 REV B

TRANSCEND® 365 miniCPAP™

Transcend 365miniCPAP kit


- Transcend 365 miniCPAP™
- Secure padded travel bag with storage compartments
- AC power adapter with international power cords
- USB cord
- Quick start guide



new TRANSCEND® 365 miniCPAP™



the
**SMART
CPAP**
at home & on the go

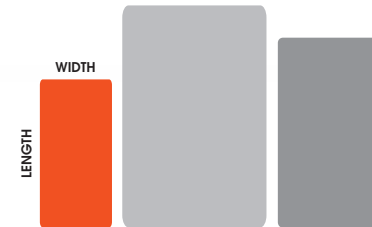

The
Smart CPAP
at home and
on the go.


Footprint size comparison

	365 miniCPAP™	DreamStation™	AirSense 10™
CPAP + HH dim. (cm)	19.9 x 9.7 x 11.5	29.7 x 19.3 x 8.4	25.4 x 15 x 11.7
CPAP + HH size (cm ³)	2,219.8	4,815	4,457.7
CPAP + HH Weight (g)	770	1,996	1,247
CFV™ Humidification	Yes	No	No
Portable Battery Option ²	Yes	No	No
Wireless Compliance Monitoring	Yes	Yes	Yes
Fully Featured	Yes	Yes	Yes

DreamStation™ is a trademark of Philips Respironics®, AirSense 10™ is a trademark of ResMed®, © 2019 Somnetics International, Inc. Protected by U.S. and foreign patents issued and pending. Transcend® 365 miniCPAP™ is a registered trademark of Somnetics International, Inc.

new TRANSCEND® 365 miniCPAP™

the
**SMART
CPAP**
at home & on the go

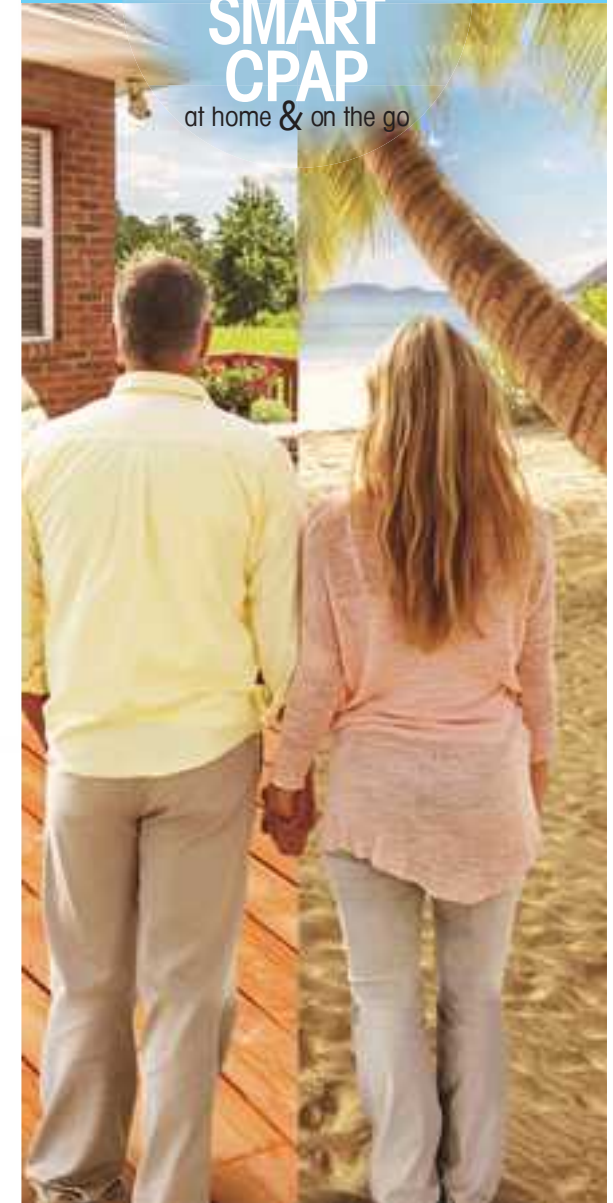

www.MyTranscend.com

WORLD'S LIGHTEST INTEGRATED CPAP + HEATED HUMIDIFIER.



New Color LCD Screen

You can now view your compliance data and settings on your 365 miniCPAP device.



Large, Tactile Buttons

Large buttons make it easy for users to navigate the menu screens and begin therapy.

The Only CPAP + Humidifier of its Kind

The 365 miniCPAP is the first and only CPAP system to use patented Capillary Force Vaporizer™ (CFV) technology. CFV delivers warm, moist air to you only when you inhale to provide comfortable therapy when you need it. The 365 miniCPAP uses about half the water of a standard heated humidifier and much less energy. This small, lightweight package can deliver full humidification almost immediately and can even operate on battery power!

NEW Battery-Powered CPAP and Humidifier

Experience the comfort of heated humidification every night, everywhere. Thanks to smart, patented power management technology, the small but powerful P10 CPAP battery can power the Transcend 365 miniCPAP AND heated humidifier so you can sleep comfortably every night, everywhere.⁽²⁾ The 365 miniCPAP conveniently delivers comfortable CPAP therapy at home and on the go in a sleek, modular design. That's smart.



Perfectly Stacked

Thoughtful design allows you to neatly stack the battery under the CPAP and humidifier for a minimal profile on the nightstand.

(2) Battery operation time and humidification level will vary depending upon therapy pressure, breath rate, tidal volume, leak and altitude.



Full-Size Features in a Travel-Size CPAP

Transcend 365 miniCPAP™ takes portable CPAP to a new dimension by introducing an incredibly small, truly portable, and unique heated humidifier. Now you can enjoy all the features and benefits of the CPAP you've come to rely on at home in a package that's half the size⁽¹⁾ and is perfect for traveling.

(1) The Transcend 365 miniCPAP is less than half the physical dimensions of the AirSense 10™ and DreamStation™. Capillary Force Vaporizer™ (CFV) is a trademark of Vapores, LLC and holds the following patents: 6634864, 7431570, 7920777, 8201752, and 792644.



new
TRANSCEND® 365 miniCPAP™
www.MyTranscend.com